



Flat  
Knit  
THERAPY



**VENOSAN®**

Compression Stockings

# VENOSAN® flat knit

## Measuring Instructions - Lower Extremity

### NOTE:

- All measurements in cm
- **c** = circumference, **l** = length
- Unless specified, all length measurements should be longitudinal and should not follow body contours

### Leg Measurements

Measuring points B-D can be done with patient sitting or lying.  
Above E to be measured with patient standing.

- cG** Two fingers/4cm below Gluteal Fold for circumference
- lG** Two fingers/4cm below Gluteal Fold to ground for length
- cF** Mid Thigh for circumference
- lF** Mid Thigh to ground for length
- cE** Top of kneecap for circumference to avoid being tight at knee pit area
- lE** Middle of Kneecap to ground for length
- cD** Two fingers/4cm below Popliteal Fossa for circumference
- lD** Two fingers/4cm below Popliteal Fossa to ground for length
- cC** Widest part of Calf for circumference
- lC** Widest part of Calf to ground for length
- cB1** Where Calf starts for circumference  
- approximately 6-8cm above point B
- lB1** Where Calf starts to ground for length
- cB** Smallest part of the ankle, above Malleolus for circumference
- lB** Smallest part of the ankle, above Malleolus to ground for length
- cY** Around Instep and Heel with foot dorsiflexed for circumference

### Foot Measurements

- cA** Around base of toes for circumference.  
Always in standing position, weight bearing.
- lA** **Horizontal finish:** Back of heel to base of toes for desired length  
**Oblique finish:**
  - lA<sup>medial</sup>** From the back of the heel to the base of large toe
  - lA<sup>lateral</sup>** From the back of the heel to the base of small toe
- lZ** **Closed Toe:** Back of heel to tip of longest toe  
(Total length of foot)

### Panty Section Measurements

Ask patient to put their hands on their waistline

- cT** Waist circumference
- lT** Waist to ground for length
- cH** Largest circumference around the hips or stomach

Measure the following length along contours of body

- lGT** Waist to 5cm/three fingers below Gluteal Fold
- lK<sup>1T</sup>** Front length from centre of waist (T) to 5cm above Gluteal Fold applying gentle pressure
- lK<sup>2T</sup>** Length at the back from the waist (T) over the buttocks to the Gluteal Fold while applying gentle pressure and following contours

### Options

**Silicone band:** Helps keep garment in place.

**Elastic band:** makes top edge more firm. Does not prevent garment from slipping down. No length adjustments required for silicone or elastic bands.

**T-Heel:** can prevent rubbing or a crease indentation on the anterior ankle. Standard option in VENOSAN Strong garments.

**Oblique G:** Garment is finished on an angle.

**Ankle inserts:** silicone ankle inserts can be requested. No additional measurements needed.

**Curved knee:** Knee area knitted in a bent position.

**Stocking in 2 parts:** No additional measurements needed. The overlap will provide no additional pressure. Overlap is 7cm in length, if you do not state otherwise.



Flat  
Knit  
THERAPY

Send completed form to Medical Accessories of Australia:  
Email: sales@medicalaccessories.com.au  
Fax: (07) 3870 5944  
Phone: 1300 527 127




VENOSAN® SOFT  
VENOSAN® PRO  
VENOSAN® STRONG

## Lower Extremity Order Form

Date:

Therapist name:	Therapist phone:	Therapist email:
Patient name:	Patient phone:	Patient schema no:
Delivery details:		

### Measurements

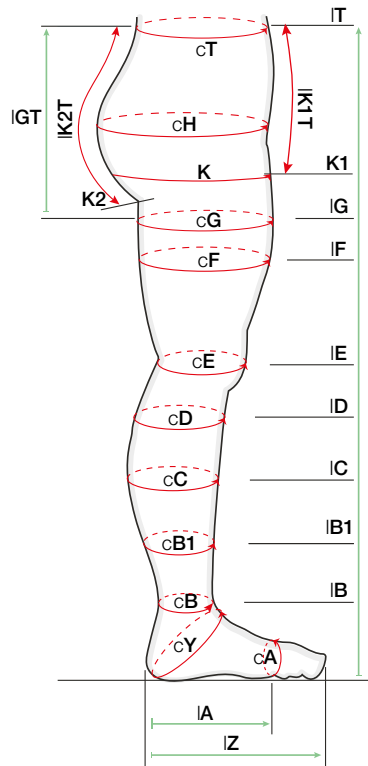


### Hip/Waist

Circumference (cm)	Length (cm)
cT	IGT
cH	IK1T
	IK2T

Foot/Leg	Left circumference (cm)	Length (cm)	Right circumference (cm)
cG	IG		cG
cF	IF		cF
cE	IE		cE
cD	ID		cD
cC	IC		cC
cB1	IB1		cB1
cB	IB		cB
cY			cY
cA	IA		cA
		IZ	
	IA-medial		} must be filled in for oblique toe
	IA-lateral		

Office use only:

Garment types:

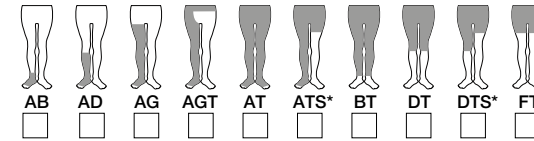

\*short leg can finish at F or G  
Additional garment styles are available. Contact Medical Accessories for details.

Product	Ccl.	Garment type	Qty. (pcs.)	Colour	
				beige	black
VENOSAN Soft	1				
VENOSAN Soft	2				
VENOSAN Pro	2				
VENOSAN Soft	3				
VENOSAN Pro	3				
VENOSAN Strong	3				
VENOSAN Strong	4				
VENOSAN Strong	4S				

### Options

#### Bands

- 3 cm silicone band
- 5 cm silicone band
- 5 cm elastic band

#### Other

- open toe
- closed toe (please fill in IA and IZ)
- oblique toe (please fill in IA-medial and IA-lateral)
- T-heel
- horizontal fly
- open crotch

- full compression body part
- oblique G \*\*
- curved knee
- stocking in 2 parts (ie. AD + DG)

\*\* With AG and DG only possible in combination with 5 cm silicone or elastic band

Comments: